## ENG (ENGLISH) - INDICATIONS

The 50 cm Disposable Radiopaque Endomyocardial Biopsy Forceps are designed for Right Ventricular Biopsies using the Jugular Approach, and for Right or Left Ventricular Pediatric Biopsies using the Femoral Approach.

The 105 cm Disposable Radiopaque Endomyocardial Biopsy Forceps are designed for for Right or Left Ventricular Biopsies using the Femoral Approach

### PRECAUTIONS

- For single use only. Do not reuse or resterilize. Sterilized by Ethylene Oxide (EiO).
- Federal (USA) law restricts this device to sale by or on the order of a physician.
- The Forcep should be thoroughly rinsed with Heparinized Saline prior to each specimen collection during the procedure.
- The Heart should be routinely monitored by EKG during the procedure.

### WARNINGS

- Inspect the packaging and product for damage prior to use. Do not use if the package is damaged as sterility may be compromised.
- The Biopsy Forcep is designed and intended for single, one patient use. DO NOT REUSE. Structural or functional integrity of the device may be compromised through cleaning, reesterilization, and subsequent reuse.

### SUGGESTED DIRECTIONS FOR USE

The Biopsy Forcep should be advanced through a Sheath Introducer or Guiding Catheter and into the Right or Left Ventricle using Fluoroscopic or Echocardiographic Guidance,

NOTE: Confirm the position of the Long Sheath Introducer or Guiding Catheter in Ventricle prior to each and after each specimen collection.

- The jaws of the Biopsy Forcep should be opened as the Forcep emerges from the distal end of the Long Sheath Introducer or Guiding Catheter.
- The open jaws of the Biopsy Forcep should be positioned at the Heart Wall, closed firmly and sufficient pressure should be maintained on the handle to assure retention of the specimen during withdrawal through the Sheath Introducer or Guiding Catheter.
- Continuous flushing of the Introducer Sheath or Guiding Catheter with Heparinized Saline should be conducted while slowly withdrawing the Biopsy Forcep/Specimen.

NOTE: Biopsy Forceps labeled as formable can be shaped at the distal tip. To form or modify the angle of the curve, grasp the distal tip between the thumb and index finger. Slide fingers along tip body while bending to the desired angle. Do not kink or form the Biopsy Forcep to more than 90 degrees.

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
## Endomyocardial Biopsy Forcep

## Multi-language Instructions For Use (IFU)

Languages	Page	Languages	Page
English	1	Hungarian	4
Russian	1	Lithuanian	5
Czech	1	Latvian	5
Danish	1	Maltese	5
Turkish	2	Norwegian	5
German	2	Dutch	6
Greek	2	Polish	6
Spanish	3	Portuguese	6
Estonian	3	Romanian	7
Finnish	3	Slovak	7
French	3	Swedish	7
Italian	4	繁體中文	7
Croatian	4		

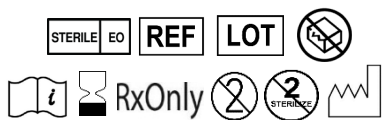**Argon Medical Devices, Inc.**  
1445 Flat Creek Road,  
Athens, Texas 75751, USA  
Tel: +1 (903) 675 9321  
Tel: +1 (800) 927 4669  
www.argonmedical.com

**EMERGO EUROPE**  
Prinsessegracht 20 2514  
AP The Hague The  
Netherlands  
+31 70 345 8570

The symbols glossary is located electronically at [www.argonmedical.com/symbols](http://www.argonmedical.com/symbols)

Made in USA  
65-2000-01/0819T