ARGON MEDICAL DEVICES, INC.
1445 Flat Creek Road, Athens, Texas 75751 USA
Tel: +1 800-927-4669;
Tel: +1 903-675-9321
www.argonmedical.com

EMERGO EUROPE
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands
+31 70 345 8570

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
Emergo Consulting (UK) Limited
c/o Cr 360 - UL International
Compass House, Vision Park Histon
Cambridge CB24 9BZ
United Kingdom

X9585357-0321

Date of Issue : 08/14/2023 Rev. B

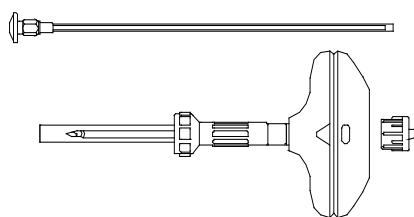MD STERILE EO RxOnly

<https://www.argonmedical.com/resources/product-information>

English	The symbols glossary is located electronically at www.argonmedical.com/symbols
Bulgarian	Речникът на символите е наличен в електронен вид на адрес www.argonmedical.com/symbols
Czech	Slovník symbolů v elektronické formě je umístěn na www.argonmedical.com/symbols
Danish	Symbolforklaringen findes elektronisk på www.argonmedical.com/symbols
German	Das Glossar der Symbole kann unter www.argonmedical.com/symbols abgerufen werden
Greek	Το γλωσσάρι των συμβόλων παρέχεται ηλεκτρονικά στη διεύθυνση www.argonmedical.com/symbols
Spanish	El glosario de símbolos puede consultarse en línea en la dirección www.argonmedical.com/symbols
Estonian	Sümbolite seletuste kogu on elektrooniliselt saadaval aadressil www.argonmedical.com/symbols
Finnish	Elektroninen symbolisanasto on osoitteessa www.argonmedical.com/symbols
French	Le glossaire des symboles est disponible au format électronique à l'adresse www.argonmedical.com/symbols
Croatian	Pojmovnik simbola nalazi se u elektroničkom obliku na www.argonmedical.com/symbols
Hungarian	A szimbólumok összefoglalása elektronikusan itt található: www.argonmedical.com/symbols
Italian	Il glossario dei simboli è disponibile in formato elettronico all'indirizzo www.argonmedical.com/symbols
Latvian	Simbolių žodyno elektroninė forma pateikta www.argonmedical.com/symbols
Lithuanian	Simbolų glosarijs ir atrodams elektroniskā veidā vietnē www.argonmedical.com/symbols
Dutch	Uitleg over de symbolen is elektronisch beschikbaar op www.argonmedical.com/symbols
Norwegian	Symbolforklaringen befinner seg elektronisk på www.argonmedical.com/symbols
Polish	Glosariusz symboli w postaci elektronicznej dostępny jest pod adresem www.argonmedical.com/symbols
Portuguese	O glossário de símbolos está localizado eletronicamente no endereço www.argonmedical.com/symbols

Romanian	Glosarul simbolurilor este stocat electronic la www.argonmedical.com/symbols
Russian	Электронный глоссарий символов расположен по адресу www.argonmedical.com/symbols
Slovak	Glosár symbolov je v elektronickej podobe k dispozícii na stránkach www.argonmedical.com/symbols
Slovenian	Slovar simbolov se nahaja v elektronski obliki na www.argonmedical.com/symbols
Swedish	Symbolforklaringen hittas elektroniskt på www.argonmedical.com/symbols
Turkish	Semboller sözlüğü elektronik ortamda www.argonmedical.com/symbols adresinde yer almaktadır

ENGLISH

Osty-Core™ Bone Biopsy Needle

Intended Use/Purpose: The Osty-Core™ Bone Biopsy Needle is used to harvest bone marrow specimens.

Device Description: Osty-Core™ Bone Biopsy Needle is designed with unique tip geometry to provide superior needle entry and coring characteristics without distorting the biopsy specimen. The ergonomic handle features a twist-lock design that stabilizes the needle during placement, while a Luer lock connector on the handle provides secure syringe attachment. It has an adjustable, removable needle stop to safely control the depth of penetration and can be used with co-axial introducer needles.

Indications For Use: Intended for the purpose of obtaining a core biopsy of bone.

Contraindications: For use only for obtaining a core biopsy of bone as determined by a licensed physician. Only a physician familiar with the possible side effects, typical findings, limitations, indications and contraindications of bone biopsy should use these needles. Physician judgement is required when considering biopsy on patients with bleeding disorder or receiving anti-coagulant medications.

Potential Complications/Risks: It is not uncommon for patients to experience the following :

- Nerve injury including puncture of the spinal cord or nerve roots potentially resulting in radiculopathy, paresis or paralysis.
- Thrombus or other materials resulting in symptomatic pulmonary embolism or other clinical sequelae.
- Hemothorax or pneumothorax.
- Infection.
- Unintended puncture wounds including vascular puncture and dural tear.
- Bleeding or hemorrhage
- Hematoma
- Pain

When used as intended by a physician, the benefits of using the Bone Marrow Biopsy Needle outweigh the risks associated with the use of the device.

*** Pulmonary and cardiovascular related complications occur mostly in rare situations when the device is used for sternal aspirations.

Warning

- Inspect the package integrity before use.
- Do not use if package is open or damaged and if the expir date has been exceeded.
- Do not continue to use if any of the component are damaged during procedure.
- Do not bend the needle excessively, it may result damage/breakage to the device.
- This device was designed, tested and manufactured for single patient use only.
- When used repeatedly in the same patient, the device should be inspected for damage or wear after each core sample is taken. Reuse or reprocessing has not been evaluated and may lead to its failure and subsequent patient illness and/or injury. DO NOT reuse, reprocess or re-sterilize this device.
- Reconditioning, refurbishing, repair, modification, of the device to enable further use is expressly prohibited.

Precaution

- This device should be manipulated only while under fluoroscopic observation with radiographic equipment that provides high quality images.
- Access to the vertebral body via pedicle requires a minimum pedicle width of 5mm.
- Knowledge of the dimensions of the site of insertion as assessed by MRI, CT, or other imaging method.

NOTE: These instructions are not intended to define or suggest any medical or surgical technique. The individual physician is responsible for the proper procedure and techniques to be employed with this device.

RECOMMENDED TECHNIQUE FOR SINGLE STICK BIOPSY:

1. Using sterile technique, prepare the skin with antiseptic and drape.
2. After having determined the most suitable point of entry, the entire area between the skin and periosteum is anesthetized. The periosteum must be well anesthetized also. The fine anesthesia needle remains in position as a guide for a parallel puncture with the Osty-Core™ Bone Biopsy Needle. The position of the needle can be controlled using Fluoroscopic or "CT" guidance.
3. Hold the Osty-Core™ Bone Biopsy Needle with the proximal end in palm and the index finger against the shaft near the tip. This position stabilizes the needle and allows better control.
4. Introduce the Osty-Core™ Bone Biopsy Needle through the skin and bring it into contact with the area of biopsy interest.
5. The fine anesthesia needle can now be removed.
6. Using gentle, but firm pressure, advance the needle, rotating it in an alternating clockwise/counterclockwise motion through the cortex/area of interest. Once the cannula passes through the cortex, the stylet may be removed.
7. After reaching the lesion/area of interest, attach a syringe with luer taper to the hub of the Osty-Core™ Bone Biopsy Needle using a firm push/twist motion.
8. Apply suction by withdrawing the syringe plunger. With vacuum applied, slowly rotate and remove the needle cannula with specimen.
9. Remove the syringe from the needle cannula.
10. Push the specimen out of the cannula using the obturator.

RECOMMENDED TECHNIQUE FOR MULTIPLE BIOPSY SPECIMENS--PROCEDURE FOR CO-AXIAL BONE BIOPSY:

- Using sterile technique, prepare the skin with antiseptic and drape.
- After having determined the most suitable point of entry, the entire area between the skin and periosteum is anesthetized. The periosteum must be well anesthetized also. The fine anesthesia needle remains in position as a guide for a parallel puncture with the Osty•Core™ Bone Biopsy Needle. The position of the needle can be controlled using Fluoroscopic or "CT" guidance.
- A 15ga. Co-Axial Introducer Needle (refer to accompanying table) is placed through the surrounding tissue to the edge of the bone lesion/area of interest. Once in position, the stylet should be removed.

Osty•Core™ Bone Biopsy Needle	Co-Axial Introducer Needle
OBN16075	N/A
OBN16125	MCXS1610AX
OBN16165	MCXS1616AX

- Hold the Osty•Core™ Bone Biopsy Needle with the proximal end in palm and the index finger against the shaft near the tip. This position stabilizes the needle and allows better control.
- The Osty•Core™ Bone Biopsy Needle is then placed through the co-axial introducer needle and into contact with the area of biopsy interest.
- The fine anesthesia needle can now be removed.
- Using gentle, but firm pressure, advance the Osty•Core™ Bone Biopsy Needle, rotating it in an alternating clockwise/counterclockwise motion through the cortex/area of interest. Once the cannula passes through the cortex, the stylet may be removed.
- After reaching the lesion/area of interest, attach a syringe with luer taper to the hub of the Osty•Core™ Bone Biopsy Needle using a firm push/twist motion.
- Apply suction by withdrawing the syringe plunger. With vacuum applied, slowly rotate and remove the needle cannula with specimen.
- Remove the syringe from the needle cannula.
- Push the specimen out of the cannula using the obturator.
- For additional specimens, repeat as required from step #7.
- After obtaining desired specimen(s), remove the co-axial introducer needle cannula.

Disposal: After use, this product may be a potential biohazard. Handle in a manner which will prevent accidental puncture. Dispose in accordance with applicable laws and regulations.

Storage: Store at controlled room temperature.

NOTE: In the event a serious incident related to this device occurs, the event should be reported to Argon Medical at quality.regulatory@argonmedical.com as well as to the competent health authority where the user/patient resides.

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