

CLEANER Vac® Thrombectomy System in 200 Patients

88%

Mean clot
clearance

94%

Completed in a
single session

92%

Achieved 75% or greater
thrombus resolution

0

Major complications
reported

68 min

Average procedure duration

8 min

Mean device utilization time

2.4

Device passes per patient

2.2

Venous segments treated

Study Overview

- 200 patients, retrospective
- 26 U.S. institutions
- 64.5 years mean patient age
- 52% / 48% female / male
- 65% / 35% acute / subacute thrombus
- VTE symptomatic venous thromboembolism




CLEANER Vac Thrombectomy System

PROCEDURE DETAIL & CLINICAL INSIGHT

How the Procedures Were Performed

CLEANER Vac was used as the **primary aspiration device** in patients presenting with symptomatic venous thromboembolism (VTE) across multiple participating centers. Device activation and aspiration were performed under fluoroscopic guidance until thrombus clearance was achieved, with adjunctive therapies (such as balloon angioplasty or thrombolytics) used as needed at the operator's discretion. Technical success, procedure duration, and peri-procedural safety were recorded for all cases.

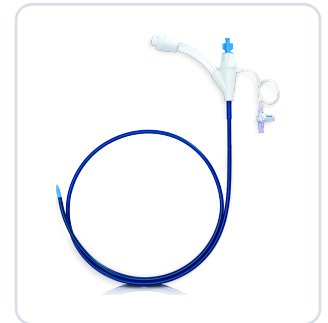

Clinical Insight

Sites with greater procedural experience demonstrated **substantially reduced blood loss**, suggesting opportunities for further optimization through technique standardization and physician training.

Ordering Information

CATALOG #	CLEANER VAC SYSTEM (BOX OF 1)
70027518	CLEANER Vac Thrombectomy System

REFERENCE

Stibbs, P. *Retrospective Multicenter Experience with the CLEANER Vac Thrombectomy System: An Analysis of 200 Patients* [abstract]. Accepted for presentation and for publication in CardioVascular and Interventional Radiology, 2026.

Abstract accepted for presentation at the CIRSE 2026 meeting; full manuscript not yet published.

RX Only.

Indication for Use: The CLEANER Vac Thrombectomy System is indicated for the removal of fresh, soft thrombi and emboli from the vessels of the peripheral venous vasculature, and for the infusion of physician-specified fluids, including thrombolytics. The CLEANER Vac Thrombectomy System is not intended for use in the pulmonary vasculature for treating of pulmonary embolism.

Contraindications for Use: Not intended for use in the coronaries and neurovasculature.

Refer to the Instructions for Use for complete warnings, and precautions.

© 2026 Argon Medical Devices, Inc. All rights reserved. Argon logo, Argon Medical Devices, and Cleaner Vac are registered trademarks of Argon Medical Devices, Inc.