

# Complete Drainage System

**SKATER™**



**The SKATER™ Drainage System** provides a complete portfolio of catheters and accessories for all-purpose, nephrostomy, biliary and centesis procedures.

### All-Purpose Drainage Catheters

For drainage of abscesses and fluid collections

- Choice Lock™ trocar stylet allows for single-handed abscess drainage placement
- Available in locking and non-locking configurations
- Expansive size offering of 6F - 16F
- Can be used safely with ethyl alcohol and doxycycline

### Kit Components Contain One:

- SKATER catheter
- Removable pigtail straightener
- Metal stiffening cannula
- Flexible stiffening cannula (Not available on 6F and 7F)
- Locking system with snap-off tab (Not available on non-locking pigtail catheters)

### Nephrostomy Drainage Catheters


For drainage of the urinary tract

- Designed for Seldinger technique
- Expansive size offering of 6F - 14F
- Can be used safely with ethyl alcohol and doxycycline

### Biliary Drainage Catheters

For drainage of the biliary tract

- Radiopaque marker band provides clear view of the most proximal drainage hole for accurate placement in the biliary tract
- Additional drainage holes for biliary duct to duodenum drainage
- Can be used safely with ethyl alcohol and doxycycline

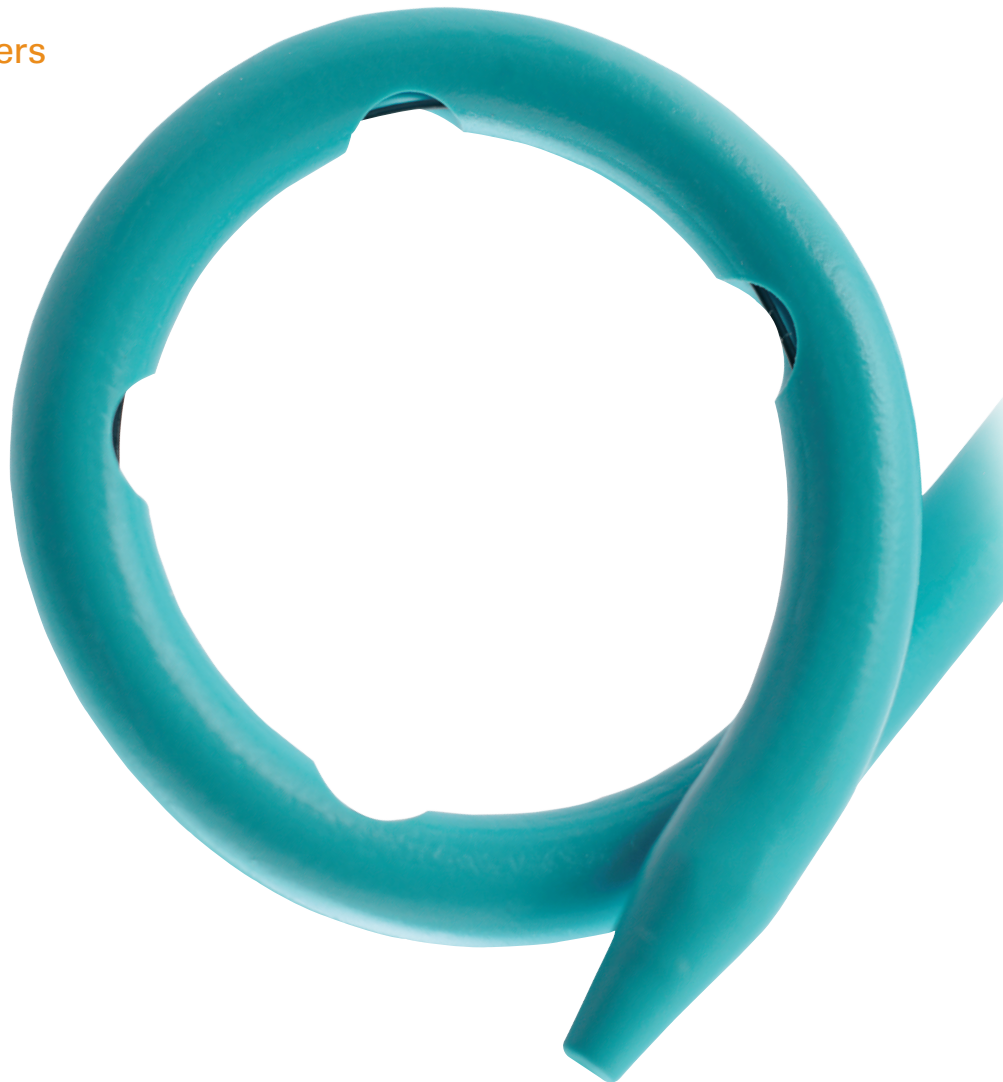

# SKATER™

# Formulated to promote

**Large skived holes**  
help maximize drainage

**Position markers**  
confirm catheter  
position after  
placement


**Locking  
shoulder tip**  
reduces buckling  
during insertion

**Removable pigtail**

**SLIP-COAT™ hydrophilic coating**  
significantly reduces friction for  
easier catheter insertion



# maximum drainage flow

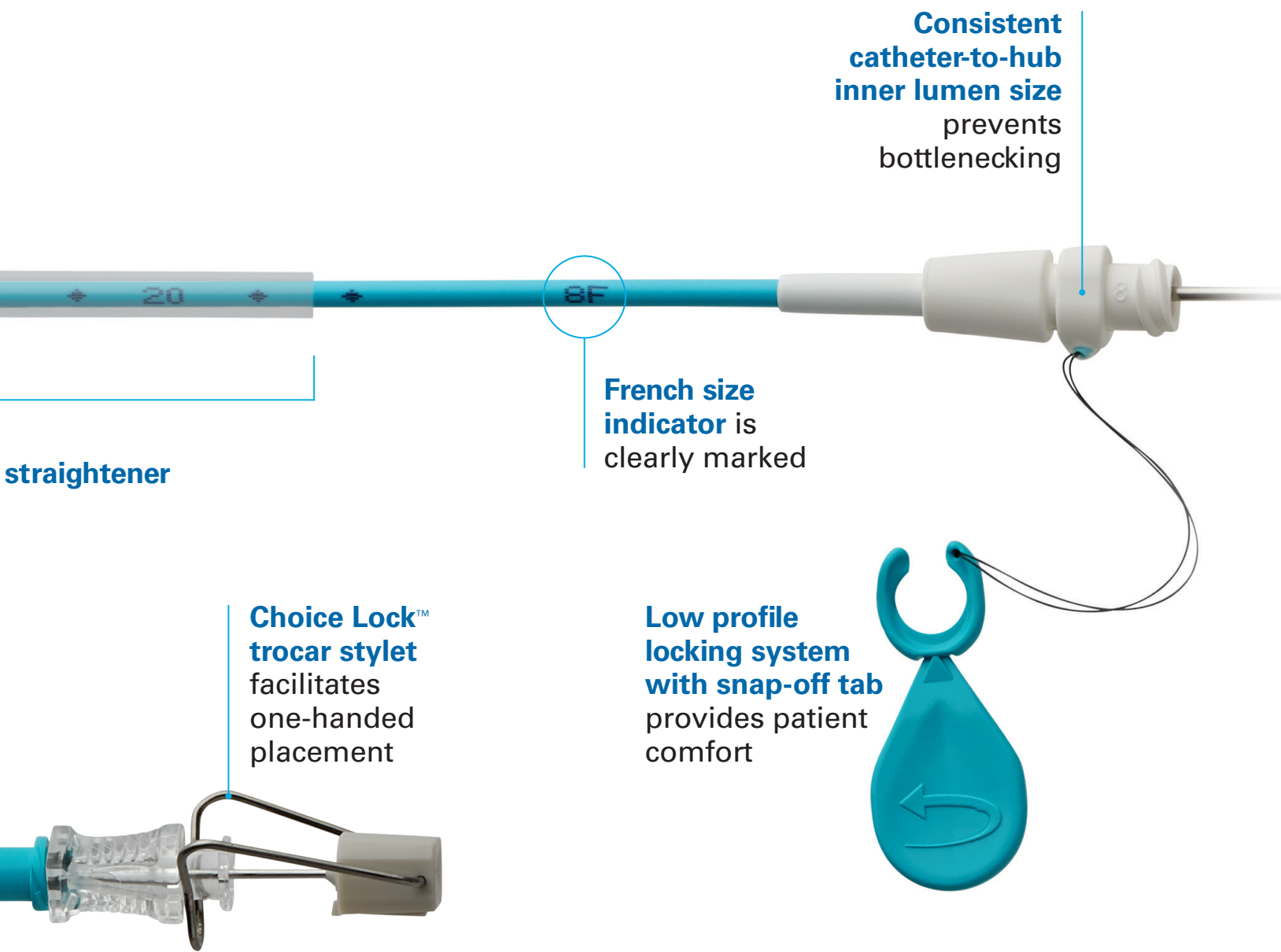

**Consistent catheter-to-hub inner lumen size** prevents bottlenecking

**French size indicator** is clearly marked

**Choice Lock™ trocar stylet** facilitates one-handed placement

**Low profile locking system with snap-off tab** provides patient comfort

straightener

## SKATER Accessories

### Centesis Catheters

- Needle with echogenic tip for precise placement under imaging
- Tapered catheter tip ensures smooth insertion
- Four radial side holes facilitate drainage within small cavities
- Hollow stylet needle permits aspiration of fluid prior to needle removal
- Vent plug allows air to exit catheter during insertion and protects against fluid leakage




### SKATER Introducer Systems

- Echogenic needle tip allows for accurate placement with ultrasound guidance
- Embedded radiopaque marker enhances visualization for accurate placement
- 4F dilator/6F sheath assembly with locking stiffening cannula facilitates easy over-the-wire placement
- Large end hole design of the sheath accepts up to both .018" and .038" guidewires simultaneously
- Available with a stainless steel or nitinol .018" guidewire



### Hawkins™ Blunt Needle Access Systems

- The adjustable needle stop localizes the outer cannula to the target, preventing forward migration
- The outer cannula has an echogenic tip, promoting accurate placement under ultrasound guidance
- Centimeter markings provide reference for accurate depth placement
- Accommodates most guidewires

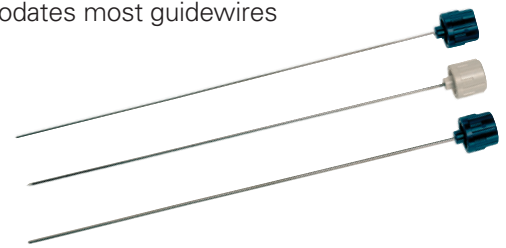

### Drainage Bags

- Kink-resistant tube can be cut and fit to the appropriate length
- Flexible tubing allows for patient mobility during drainage
- Twist-type valve provides quick and easy emptying of drainage bag
- Anti-reflux valve prevents retrograde contamination
- Fully adjustable soft belt provides patient convenience



### Introducer Sheaths/Needles (ISN)

- 5F radiopaque sheaths for guidewire insertion or fluid aspiration
- 19G needle with trocar-style tip for easy insertion
- Available in 15 cm and 20 cm lengths

## Ordering Information

SKATER All-Purpose Drainage Catheters w/ Locking Pigtail (Box of 5)	
CATALOG #	DESCRIPTION
756506020*	6F x 20 cm, accepts up to .035" guidewire
756507020*	7F x 20 cm, accepts up to .035" guidewire
756508025	8F x 25 cm, accepts up to .038" guidewire
756508030	8F x 30 cm, accepts up to .038" guidewire
756510025	10F x 25 cm, accepts up to .038" guidewire
756510030	10F x 30 cm, accepts up to .038" guidewire
756512025	12F x 25 cm, accepts up to .038" guidewire
756514025	14F x 25 cm, accepts up to .038" guidewire
756516025	16F x 25 cm, accepts up to .038" guidewire

\*Flexible stiffening cannula not available on 6F and 7F.

SKATER All-Purpose Drainage Catheters w/ Non-Locking Pigtail (Box of 5)	
CATALOG #	DESCRIPTION
756006020*	6F x 20 cm, accepts up to .035" guidewire
756007020*	7F x 20 cm, accepts up to .035" guidewire

\*Flexible stiffening cannula not available on 6F and 7F.

SKATER Nephrostomy Drainage Catheters w/ Locking Pigtail (Box of 5)	
CATALOG #	DESCRIPTION
755606025*	6F x 25 cm, accepts up to .035" guidewire
755608035	8F x 35 cm, accepts up to .038" guidewire
755610035	10F x 35 cm, accepts up to .038" guidewire
755612035	12F x 35 cm, accepts up to .038" guidewire
755614035	14F x 35 cm, accepts up to .038" guidewire

\*Flexible stiffening cannula not available on 6F.

SKATER Biliary Drainage Catheter w/ Locking Pigtail (Box of 5)	
CATALOG #	DESCRIPTION
755308040	8F x 40 cm, accepts up to .038" guidewire
755310040	10F x 40 cm, accepts up to .038" guidewire
755312040	12F x 40 cm, accepts up to .038" guidewire

SKATER Centesis Catheters (Box of 5)	
CATALOG #	DESCRIPTION
700304007	4F x 7 cm
700304010	4F x 10 cm
700304015	4F x 15 cm
700305007	5F x 7 cm
700305010	5F x 10 cm
700305015	5F x 15 cm
700305020	5F x 20 cm

SKATER Introducer Systems (Box of 5)	
CATALOG #	DESCRIPTION
651006000	4F dilator/6F x 18 cm co-axial introducer, .018" x 80 cm stainless steel guidewire w/ platinum tip, 21G x 15 cm trocar needle w/ echogenic tip
651006300	4F dilator/6F x 18 cm co-axial introducer, .018" x 80 cm nitinol guidewire w/ platinum tip, 21G x 15 cm trocar needle w/ echogenic tip
651506000	4F dilator/6F x 20 cm co-axial introducer, .018" x 80 cm stainless steel guidewire w/ platinum tip, 21G x 15 cm trocar needle
651506300	4F dilator/6F x 20 cm co-axial introducer, .018" x 80 cm nitinol guidewire w/ platinum tip, 21G x 15 cm trocar needle

Introducer Sheaths/Needles (Box of 10)	
CATALOG #	DESCRIPTION
ISN1915	19G x 15 cm w/ 5F sheath, accepts up to .038" guidewire
ISN1920	19G x 20 cm w/ 5F sheath, accepts up to .038" guidewire

Hawkins Blunt Needle Access Systems (Box of 5)	
CATALOG #	DESCRIPTION
216057	16G x 5 cm
216107	16G x 10 cm
216157	16G x 15 cm
218057	18G x 5 cm
218107	18G x 10 cm
218157	18G x 15 cm
218207	18G x 20 cm

Drainage Bags (Box of 5 or 20)	
CATALOG #	DESCRIPTION
DBAG600	600 mL drainage, w/ adjustable drain tube length, box of 20
DBAG600H	600 mL drainage, w/ adjustable drain tube length, box of 5

© 2019 Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries. All other trademarks are the property of their respective owners.

2600 Dallas Parkway, Suite 440, Frisco, TX 75034  
P: 800.927.4669 F: 903.677.9396 argonmedical.com



96-6300-14/0619C