




Reduce risk of catheter dislocation without suturing

SKATER® FIX is a latex free catheter fixation device. This consists of a central thin carrier which has hydrocolloid adhesive and is covered with non woven medical grade adhesive. The soft foam support block and velcro prevent the catheter from kinking.

The catheter will be mechanically fixed with velcro and minimal tape is involved directly on the catheter.

The velcro cover can be opened for puncture site inspection. In this way you have the possibility to look for infection or leakage whenever required.

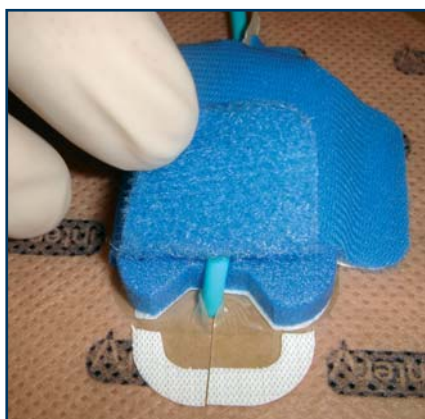

Features and benefits:

- ✓ Latex free
- ✓ Safe anchoring of drainage catheter up to one week
- ✓ Allows inspection of puncture site without change of dressing
- ✓ Easy to change
- ✓ Hydrocolloid carrier fits tightly around the catheter and will absorb body fluids from puncture site
- ✓ Minimal risk of catheter dislocation during change of fixation disk
- ✓ Foam bridge prevents catheter kinking
- ✓ Clean catheter without tape residuals
- ✓ Adjustable to catheters from 5 to 16 French
- ✓ Minimum risk of skin irritation
- ✓ Disc follows patient contour and movements



Ordering Information

SKATER® FIX (Box of 25)

Catalogue no.	Description
8025 00 001	SKATER® FIX fixation device