



The Standard in Portal Vein Access

TRAVELERTM

PORTAL VEIN ACCESS SERIES


ARGON
MEDICAL DEVICES

Your Go-To For Portal Vein Access Procedures

Go beyond the ordinary with Argon Medical's innovative **Traveler™ Portal Vein Access** Sets.

Greater Visibility.

- **Echogenic tips** enable visibility of the needle or stylet on various imaging modalities for easy detection while piercing through liver to portal vein.
- **Alignment marker** helps identify when the needle/stylet is flush with the end of the stiffening cannula.
- **Radiopaque band** on the stiffening cannula provides amplified visibility under fluoroscopy to assist with locating and orienting each component.

Improved Durability.


- **PEEK material needle/stylet catheters** help reduce the risk of fishmouthing and buckling when advancing through cirrhotic parenchyma compared to traditional PTFE.



Alignment Marker



PEEK Needle/Stylet Catheter

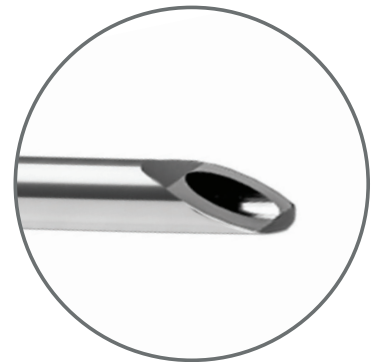

Echogenic Tips

Traveler™ Portal Vein Access Series

Argon's **Traveler Portal Vein Access Series** helps evolve portal vein access procedures through unique visibility and durability enhancements while maintaining common techniques.



TRAVELERTM16
PORTAL VEIN ACCESS SET

The Traveler 16 set consists of a blue handle, a white catheter, and a 16-gauge curved needle. The needle is shown inserted into the catheter, which is curved to follow the path of the needle.

16ga Curved Needle



TRAVELERTM21
PORTAL VEIN ACCESS SET

The Traveler 21 set consists of a blue handle, a white catheter, and a 21-gauge needle. The needle is shown inserted into the catheter, which is curved to follow the path of the needle.

21ga Needle



TRAVELERTM38
PORTAL VEIN ACCESS SET

The Traveler 38 set consists of a blue handle, a white catheter, and a 0.038-inch stylet. The stylet is shown inserted into the catheter, which is curved to follow the path of the stylet.

0.038" Stylet

Access Set Components

TRAVELER[™]16 PORTAL VEIN ACCESS SET

Catalog #TPS003

16ga x 57cm Curved Needle
7F x 51cm PEEK Needle Catheter
5F x 80cm MPA Catheter
Bending Tool
Spacer Clip



TRAVELER[™]21 PORTAL VEIN ACCESS SET

Catalog# TPS002

14ga x 52cm Stiffening Cannula
7F x 52cm Cannula Sheath
21ga x 72cm Straight Needle
5F x 67cm PEEK Needle Catheter
5F x 80cm MPA Catheter
Bending Tool
Spacer Clip



TRAVELER[™]38 PORTAL VEIN ACCESS SET

Catalog# TPS001

14ga x 52cm Stiffening Cannula
7F x 52cm Cannula Sheath
0.038" x 69cm Stylet
5F x 66cm PEEK Stylet Catheter
5F x 80cm MPA Catheter
Bending Tool
Spacer Clip



Indication

Traveler[™] Portal Vein Access Sets are intended for transjugular liver access in diagnostic and interventional procedures.

Contact your local Argon Medical representative or customer service at customer.service@argonmedical to learn more about these innovative products.

© 2022 Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries.

7800 Dallas Parkway, Suite 200, Plano, TX 75024
P: 800.927.4669 F: 903.677.9396 argonmedical.com

MADE IN THE USA

96-6500-29/0123A